

# VERTICA III VERTICAL

Architectural

## SERVICE COLUMN

arq



### -----Structure-----

#### MATERIALS

##### ● Structure:

Aluminium profile extrusion 6063 T5 alloy meeting ASTM B221-02/AA. 2 lateral internal channels duct type and front isolated open panel channel to medical gas, electrical, data and lighting. Solid state separation (not removable). Add galvanized steel pipe to wiring.

#### FINISHING

##### ● Structure:

10 microns natural anodized.

##### ● Removable Front Panels:

- 1) Aluminium sheet powder coating finishing or wood imitation by sublimated process.
- 2) MDF wood panel back support 16mm thickness with plastic laminated decorative sheet 1mm thickness. (Formica, Ralph Wilson, Nevamar covers).
- 3) MDF wood panel back support 16mm thickness with childs decorative front panel, laser color print over plastic acrylic sheet 3mm thickness, polyvinyl flexible plastic cap edges to wood back panel support.
- 4) To other materials or dimentions, consult to the factory.

#### LATERAL INSIDE STRUCTURE PATIENT LIGHTING

T8 LED tube 60cm length, 110-220v main power 60-50 hztz , 4,300°k color light (day light/white) to environment and patient reading lighting.

Plastic extrusion transparent Polycarbonate diffusers self supported to vertical bedhead unit aluminium structure.

Button type light switch to environment and patient reading lighting.



ARQ



\*\*\* Installation must be performed by medical gas certified personnel\*\*\*


GRUPO ARIGMED S. DE R.L. DE C.V. TIZAYUCA, ESTADO DE HIDALGO, MÉXICO

www.arigmed.com Arigmed Medical Gas Equipments ventasmx@arigmed.com +52 779 79 69 370



CERTIFIED



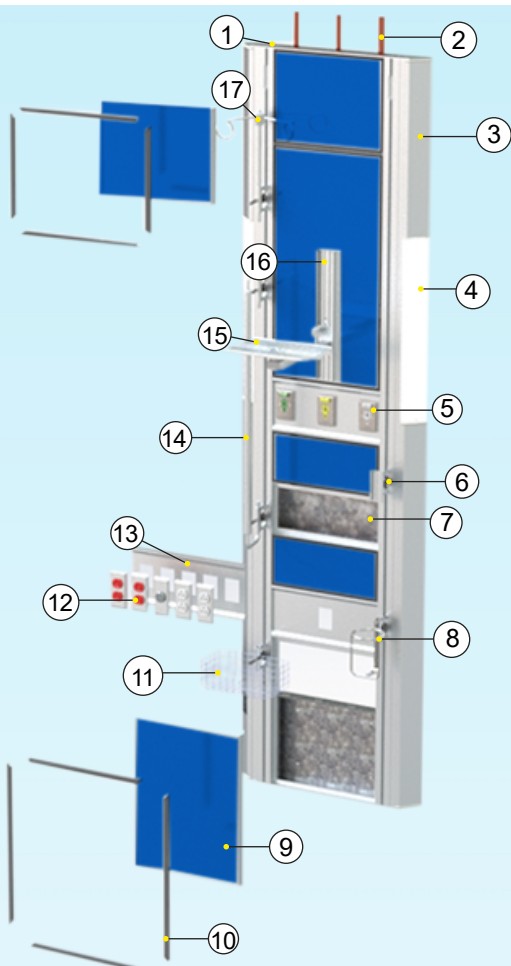

Item:	Description:
1	Aluminium top cover
2	Copper Pipe
3	Aluminium Structure frame
4	Side Patient Lighting (Inside of aluminum frame)
5	Medical Gas Outlet
6	Support and slide device to vert. rail
7	Back side cover compartment
8	1L Glass Vacuum bottle support
9	Removable front Panel
10	Plastic extrude profile screw cover
11	Medical Basket
12	Electrical Outlet
13	Aluminium service front panel
14	Infusion pump support
15	Patient monitor shelf type support
16	Vertical Service Rail
17	I.V. Support

## VERTICAL SERVICE RAILS (SUPPORT AND SLIDE)

Vertical column bedhead unit includes service rails (1 rail by each side of column) to support and slide of medical auxiliary equipment.

Complies with NFPA99-2015 exclusively with systems and devices of electrical, light and medical gases that have been assembled from factory origin.

In the case the vertical bedhead unit has been supplied partially or modified in the components and systems of factory origin, we disclaim any responsibility for risk and injuries and damages to people and properties.

## ELECTRICAL OUTLETS

It can accept all kinds of electrical outlets models according with local or international electric codes.

## MEDICAL GAS AND VACUUM INLET-OUTLETS

Mounted on the front panel, it can include a pre-assemble of rigid type K copper pipe (customer must be specify).  
vertical bedhead unit can allow the assembly of any kind of medical gas outlet.